

# Register and Save!

## Save over \$50!

Register at [Waterpik.com/Registration](https://www.waterpik.com/Registration) to receive exclusive discount code.

## We'd Love to Hear from You!

Connect with Us  
@Waterpik



Call: 1-800-525-2774  
[www.waterpik.com](http://www.waterpik.com)

## Waterpik® 14 Day Challenge App

Make water flossing a habit! Try our new mobile App for the best water flossing experience.



Offer includes  
**50% off**  
your next order  
of assorted  
tips + tip  
storage case.



waterpik® NANO PLUS  
waterflosser™

# Congratulations!

You are now on your way  
to better oral health!

### ✓ Register and Save

Register your product at [waterpik.com/registration](https://www.waterpik.com/registration) to receive over **\$50** in savings.

### ✓ Write a Review

Creating innovative Waterpik® products that inspire better oral health is our passion. We invite you to leave a review by sharing your thoughts on the Waterpik® Nano Plus Water Flosser on [waterpik.com](https://www.waterpik.com) or the retailer's website.

### ✓ We're Here to Help

We hope you are loving your new Waterpik® Nano Plus Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado. Waterpik® Nano Plus Water Flosser is backed by a 2-year manufacturer's warranty.

Call: 1-800-525-2774  
Monday-Friday  
7 am - 5 pm MST



Chat or E-mail  
[service@waterpik.com](mailto:service@waterpik.com)




**Quick Start  
Guide Inside!**




## GETTING STARTED

- 1 Fill reservoir with warm water and replace it firmly on the base.

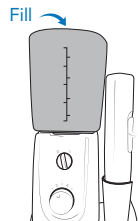


- 2 Insert tip firmly until it clicks into place.

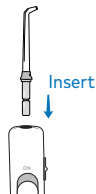

## PRIMING UNIT (First Time Use Only)


- 3 Turn pressure to high, point tip into sink, turn power ON until water flows. Turn power OFF.

Turn pressure to high, point tip into sink, turn power ON and water ON until water flows. Turn power and water OFF.

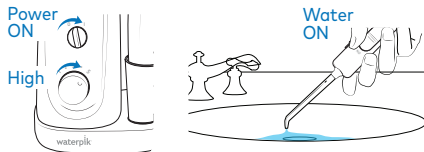


## WATER FLOSS

- 4 Place tip in mouth. Start on low pressure. Turn power ON and increase pressure to your liking.

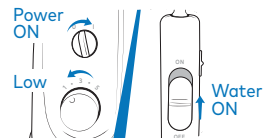


- 5 Lean over sink. Close lips enough to prevent splashing, while still allowing water to flow from mouth into sink.

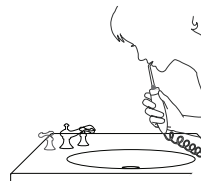

## TECHNIQUE

- 6 Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

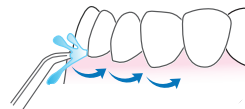

**TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.**

## NEED HELP?

Visit [Waterpik.com/WelcomeWF](https://www.waterpik.com/WelcomeWF)

To watch a How To Video and find additional product support.

**Call Us! 1-800-525-2774**