# waterpik

## Congratulations!

You are now on your way to better oral health!


**REGISTER AND SAVE**  
Register your product and receive over \$50 in savings. [waterpik.com/registration](http://waterpik.com/registration) (US) or [waterpik.ca/registration](http://waterpik.ca/registration) (CA).

**WE'RE HERE TO HELP**  
We hope you are loving your new Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado.

Call: 1-800-525-2774 (USA)  
1-888-226-3042 (Canada)  
Monday-Friday  
7 am - 5 pm MST  
Chat or E-mail  
[service@waterpik.com](mailto:service@waterpik.com)

INSTRUCTION MANUAL INSIDE →

### Waterpik® Water Flosser Models WF-16/WF-17 Series

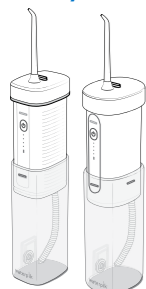**LET US HELP YOU GET STARTED!**  
We have many helpful tools to provide you with the best possible experience. Scan either QR code below or visit: [waterpik.com/slide](http://waterpik.com/slide) (US) [waterpik.ca/slide](http://waterpik.ca/slide) (Canada)

QUICK START GUIDE INSIDE →

## SAFEGUARD STATEMENTS

### READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### DANGER:

#### To reduce the risk of electrocution:

- Do not handle the user supplied USB charger with wet hands.
- Do not place or drop the user supplied USB charger into water or other liquid.
- Check the charging cord and user supplied USB charger for damage before the first use and during the life of the product.

### WARNING:

#### To reduce the risk of burns, electrocution, fire, or injury to persons:

- Use this product only as indicated in these instructions or as recommended by your dental professional.
- A user supplied USB charger with the following specifications is required to charge this product: USB Type-A, Rated Output: 5.0 VDC, 1.0 A.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord and user supplied USB charger away from heated surfaces.
- Unplug USB charger when not in use.
- Do not use the Water Flosser for more than 5 minutes in each two hour period.

- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into the nose or ear.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- Only use tips and accessories that are recommended by Water Pik, Inc.

- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with a new Pik Pocket™ tip.
- Do not drop or insert any foreign object into any opening on the product.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- There are no consumer serviceable electrical parts in this product and it does not require maintenance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

## SAVE THESE INSTRUCTIONS

SYMBOLS	
The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.	
Symbol	Definition
	Refer to instruction manual
	Manufacturer
	Ingress Protection Code
	IEC Class II Equipment
	Detachable Power Supply
	Direct Current (DC) Voltage
	Alternating Current (AC) Voltage
	Do not immerse this device in any liquid.
	Do not use this device in a shower, bathtub and do not immerse in any liquid.
	Authorized Representative in the European Community
	CE marking of conformity
	UK Conformity Assessed
	Regulatory Compliance Mark
	Separate collection for waste electrical and electronic equipment.

### Limited Two Year Warranty

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. Use with improper voltage will damage product and VOID warranty.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE

DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

### FOR USA ONLY

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

### Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

### Arbitration

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY.

Further, by purchasing, installing, or using this product you agree that arbitration must be on

an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public.

Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

### FOR CANADA ONLY

This device complies with Industry Canada's license-exempt RSSs. Operation is subject to the following two conditions: (1) This device must accept any interference, including interference that may cause undesired operation of the device.

Call: 1-800-525-2774 (USA)  
1-888-226-3042 (Canada)  
Monday-Friday  
7 am - 5 pm MST  
Chat or E-mail  
[service@waterpik.com](mailto:service@waterpik.com)

CONNECT WITH US @WATERPIK

Call: 1-800-525-2774 (US) [www.waterpik.com](http://www.waterpik.com)  
Call: 1-888-226-3042 (Canada) [www.waterpik.ca](http://www.waterpik.ca)

**waterpik**  
Water Pik, Inc., a subsidiary of Church & Dwight Co., Inc.  
1730 East Prospect Road  
Fort Collins, CO 80553-0001 USA  
1-800-525-2774  
1-888-226-3042  
[www.waterpik.com](http://www.waterpik.com)  
[www.waterpik.ca](http://www.waterpik.ca)

Designed and Tested in the USA. Made in China.


Form No. 20032124-F AA  
©2022 Water Pik, Inc.

## Quick Start Guide

### GETTING STARTED

- 1 Connect the charging cable to user supplied USB charger. Plug the charger into the wall outlet, then connect the other end of the charging cable to the unit, and charge.


**Charge Unit for 4 Hours Prior to First Use.**

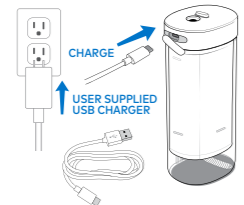
Ensure contacts on BOTH the cable and unit are clean and dry, and the reservoir is emptied and reset (close to the collapsed port) before connecting the charger.

- 2 Press tip firmly into unit until it clicks into place. To remove tip, push the tip eject button.


**Caution: Do not press tip eject button while unit is running.**

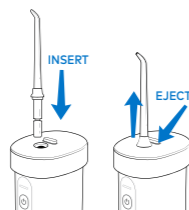
- 3 Lift the flip top on reservoir and fill with lukewarm water. Refill as needed.

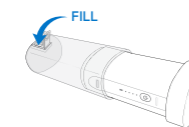- 4 Lean over sink and close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn unit ON.

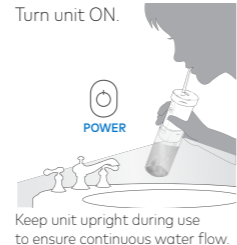- 5 Start on low pressure. Pressing the power button, while on, will change the mode. Pressing power button 4 times or holding down on power button will turn unit OFF.

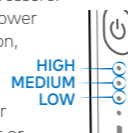- 6 Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

**TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.**

**Charging Cable**  
20030405-1 White  
WF-16/WF-17

### Replacement Tips:

#### Travel Plug

**Orthodontic Tip (OD-100E)**  
Replace every 3 months

**Plaque Seeker™ Tip (PS-100E)**  
Replace every 3 months

**Classic Jet Tip (JT-450E)**  
Replace every 6 months

Full line of replacement tips/ attachments may be purchased online at [www.waterpik.com](http://www.waterpik.com) (US) [www.waterpik.ca](http://www.waterpik.ca) (Canada), or by phone from Water Pik, Inc. at 1-800-525-2774 (US) 1-888-226-3042 (Canada).

Rated Input: 5.0 VCD, 1.0 A  
Rechargeable Battery: Lithium-ion

### Have Questions or Need More Information?

- FULL INSTRUCTION MANUAL
- CLEANING
- TIPS
- TROUBLESHOOTING
- COMMON QUESTIONS

### Let Us Help You Get Started!

We have many helpful tools to provide you with the best possible experience.

Visit: [waterpik.com/slide](http://waterpik.com/slide) (US) [waterpik.ca/slide](http://waterpik.ca/slide) (Canada)

### SERVICE MAINTENANCE

Waterpik® Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to [www.waterpik.com](http://www.waterpik.com) (US) or [www.waterpik.ca](http://www.waterpik.ca) (Canada). You may also call our toll-free customer service number at 1-800-525-2774 (US) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to [www.waterpik.com](http://www.waterpik.com) (US) or [www.waterpik.ca](http://www.waterpik.ca) (Canada)

### CLEANING

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild nonabrasive cleanser to wipe exterior. The water reservoir is removable for easy cleaning. It is recommended that the reservoir is removed and cleaned weekly.

### BATTERY RECYCLING

**Warning:** The battery of this product is not replaceable. Please do not disassemble the product unless you want to dispose it. Otherwise the warranty service will be void.

This product contains a **rechargeable Lithium-ion battery** and electronic waste. In order to protect the environment, please do not dispose along with household waste. When the battery life ends, please recycle according to local administrative regulations. Please contact local waste administrative department for any questions.

### BATTERY REMOVAL INSTRUCTIONS:

If you intend to dispose of the battery on your own, please follow the BATTERY REMOVAL INSTRUCTIONS to remove the **rechargeable Lithium-ion battery**. Contact your local waste administrative department with any questions.

- 1 The power should be off and the battery fully discharged before removal. It may take several power cycles before the battery is fully discharged.
- 2 Cut off top of unit. Cut off top of unit.
- 3 Slide battery out of chassis. Cut wires individually (avoid cutting positive and negative at the same time to prevent shorting.)

The removed battery should be disposed safely by a special recycling center.

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