


waterpik EVOLUTION™ waterflosser®


QUICK START GUIDE WF-07

Visit
waterpik.com/welcomeWF
Register your product and
receive over \$50 in savings.
waterpik.com/registration


1 First Use

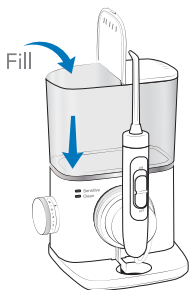

Fill reservoir with warm water and replace it firmly on the base.

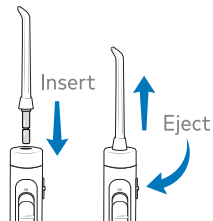

Insert tip firmly until it clicks into place.

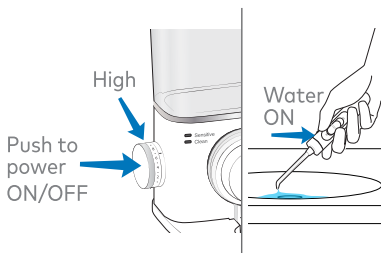


Turn pressure to high, point tip into sink, PUSH power **ON** button. Turn water switch **ON** until water flows. Turn both **OFF**.


waterpik EVOLUTION™ waterflosser®

QUICK START GUIDE WF-07


2 Water Floss

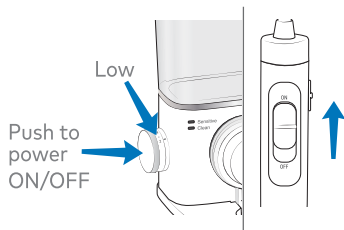

Start on low pressure. Push power **ON**. Turn water floss switch **ON** and increase pressure to your liking.

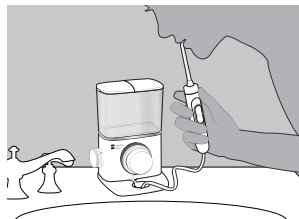

Lean over sink. Close lips enough to prevent splashing, while still allowing water to flow from mouth into sink.

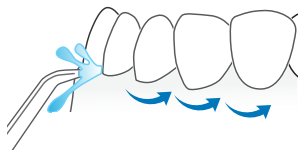

Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.