
waterpik™

Waterpik™ Complete Care Model CC-04

Soins complets Waterpik™
Modèle CC-04

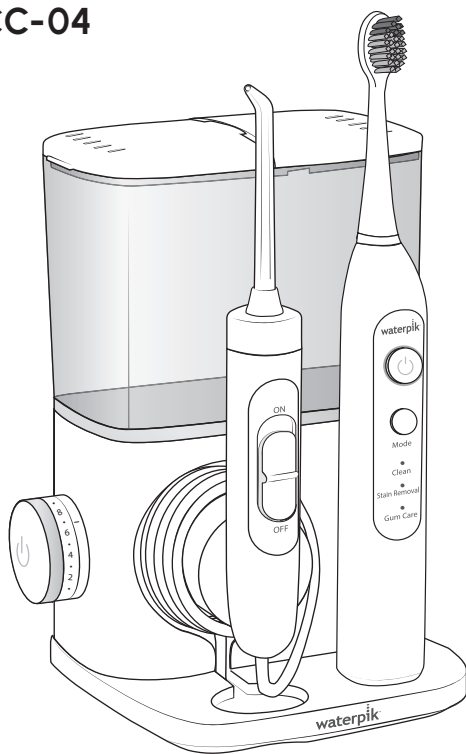www.waterpik.com

www.waterpik.ca

SAFEGUARD STATEMENTS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle plug with wet hands.
- Check the supply cord and plug for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in this Instruction Manual or as recommended by your dental professional.
- Never operate this product if it has a damaged cord or power supply, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-800-525-2774 (US) 1-888-226-3042 (Canada) or visit www.waterpik.com (US) www.waterpik.ca (Canada).
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or power supply.
- Keep electrical cord away from heated surfaces.
- Never drop or insert any object into any hose or opening.

- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Remove Water Flosser reservoir to fill. Fill reservoir with warm water or other dental professional recommended solutions only.
- Only use tips and accessories recommended by Water Pik, Inc. as identified in the instruction manual.
- The supply cord cannot be replaced. If the cord is damaged, the appliance must be scrapped.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use this product for more than five minutes in each two-hour period.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- This product does not contain any user serviceable electrical parts and it does not require maintenance.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device containing magnets or electronics, contact your physician or the device manufacturer prior to using the appliance.
- Do not expose the battery or appliance to excessive temperatures.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.

SAVE THESE INSTRUCTIONS

TABLE OF CONTENTS

Product Description	6
Getting Started	8
Care of your Waterpik™ Toothbrush	9
Cleaning and Troubleshooting your Waterpik™ Water Flosser	12
Battery Care, Battery Disposal/Recycling and Battery Removal	13-14
Limited Two-Year Warranty	15

SYMBOLS


The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.

Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Do not immerse this device in any liquid.
	Manufacturer		Do not use this device in a shower, bathtub and do not immerse in any liquid.
IPX7 or IPX4	Ingress Protection Code		Authorized Representative in the European Community
	IEC Class II Equipment		CE marking of conformity
	Detachable Power Supply		UK Conformity Assessed
	Direct Current (DC) Voltage		Regulatory Compliance Mark
	Alternating Current (AC) Voltage		Separate collection for waste electrical and electronic equipment.

PRODUCT DESCRIPTION

Water Flosser

1. Hinged reservoir lid
2. Reservoir
3. Hose
4. Water flosser pressure control dial
5. Power button (push to turn unit ON/OFF)
6. Magnetic handle cradle
7. Tip
8. Tip eject switch
9. Water ON/OFF switch
10. Charger pocket
11. Base

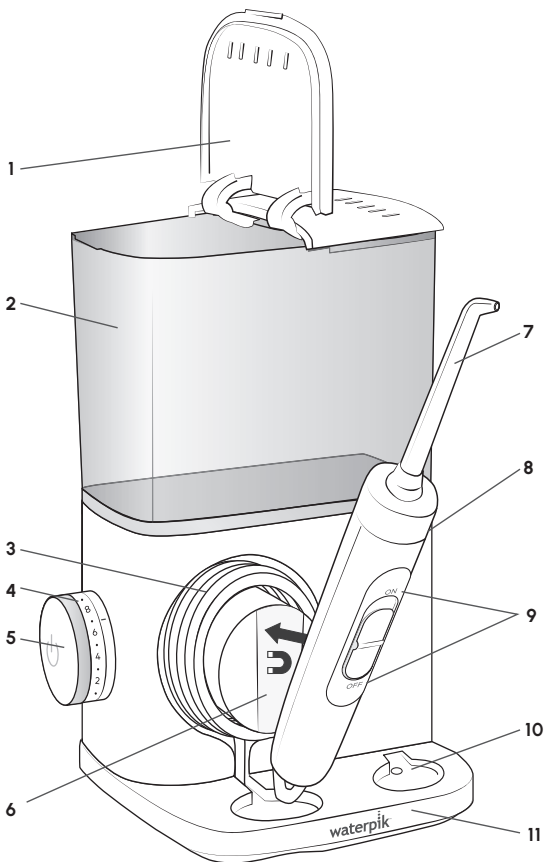TIPS

A. Orthodontic Tip*

- Braces
- General use

B. Pik Pocket™ Tip*

- Periodontal pockets
- Furcations

C. Plaque Seeker™ Tip*

- Implants
- Crowns
- Bridges
- Retainers
- General use


D. Classic Jet Tip**

- General use

E. Tongue Cleaner**

- Fresher breath

NOT ALL TIPS INCLUDED IN ALL MODELS

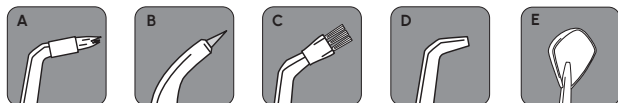*Replace Every 3 months ** Replace Every 6 months

Sensonic™ Toothbrush

12. Standard Sensonic™ Brush Head

13. LED charge indicator

14. Toothbrush power button

15. Toothbrush mode button

16. LED mode indicator

- A. Clean
- B. Stain Removal
- C. Gum Care

17. Sensonic™ Travel Case

Brush Heads and Accessories



F. Sensonic™ Toothbrush Brush Head*

- Plaque removal
- Whitening


*Replace Every 3 months

Replacement brush heads and accessories may be purchased online at www.waterpik.com (US) www.waterpik.ca (Canada), or by phone from Water Pik, Inc. at 1-800-525-2774 (US) or 1-888-226-3042 (Canada). To find detailed information about tips, other accessories, and product usage please visit our website at www.waterpik.com (US) or www.waterpik.ca (Canada).

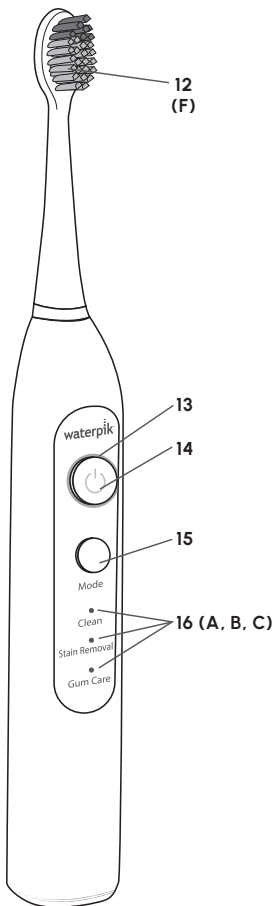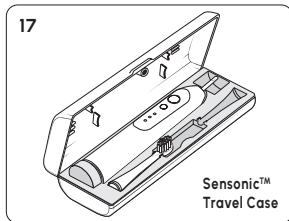Sensonic™ Toothbrush

GETTING STARTED

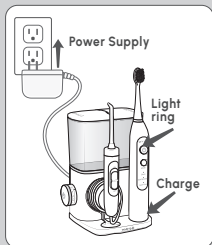**PREPARING THE UNIT****Powering the Unit**

Plug the power supply into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

Charging the Sonic™ Toothbrush

Charge for 8 hours before first use. Light ring will glow white on handle while charging. When fully charged all lights will shut off. Light ring will be solid red when battery is running low. Light ring will blink red when battery is fully depleted and needs to be recharged. A full charge will last for approximately 2 weeks of typical use (2 minutes of brushing twice a day).

Note: Charging of the toothbrush will pause temporarily when flosser is operating. Toothbrush charging will resume automatically when flosser turns off.


Modes

A. Clean: Outstanding everyday cleaning performance.

B. Stain Removal: Gently removes surface stains.



C. Gum Care: Gently enhances gum stimulation.


Note: Modes can only be changed when in operation.

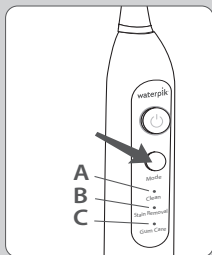**Sonic™ Brushing**

- Slide the brush head onto handle and push it down firmly.
- Place the brush head into your mouth before turning the unit on.
- Select mode by pressing mode button.
- For best results, position brush along gumline at a 45-degree angle.
- Use light pressure with a small, slow back and forth motion covering one or two teeth at a time.
- **Quad Timer:** The Toothbrush will create a short pause every 30 seconds, to remind you to move to the next quadrant of your mouth.
- **2-Minute Timer:** The brush will automatically shut off after 2 minutes of run time.

Bristles change from green to yellow when time to replace. (For optimal performance replace brush heads every 3 months.)

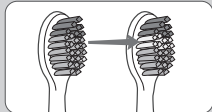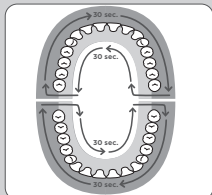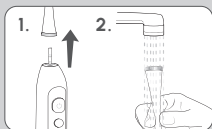Cleaning



1. **After each use**, turn off unit, remove the brush head from the handle and rinse brush head under running water.
2. Wipe the handle dry after each use.
3. Periodically, clean the charger puck, first unplugging the unit and then wiping with a damp cloth.

Service Maintenance

Waterpik™ Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to www.waterpik.com (US) or www.waterpik.ca (Canada). You may also call our toll-free customer service number at 1-800-525-2774 (US) or 1-888-226-3042 (Canada). Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Still have questions? Go to www.waterpik.com (US) or www.waterpik.ca (Canada).

Preparing the Reservoir

Remove the reservoir, lift the reservoir lid, and fill with lukewarm water. Replace the reservoir over the base and press down firmly and close the lid.

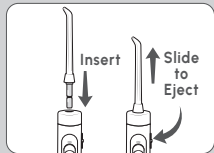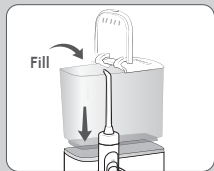Inserting and Removing Tips

Insert tip into the center of the Water Flosser handle, pressing down firmly until tip clicks into place. The colored ring will be flush with the top of the handle if the tip is correctly locked into place.

To remove tip from the handle, slide the tip eject button on the handle.



Caution: Do not use tip eject slide while unit is running.

Prime Unit (First Time Use Only)

Set pressure control dial on base to high (pressure setting No. 10), point the handle and tip into sink. Turn unit ON and water ON until water flows from the tip. Turn Power OFF.

Adjusting the Pressure Setting

Turn the pressure control dial on the base of the unit to the lowest setting (pressure setting No. 1) for first-time use. Gradually increase pressure over time to the setting you prefer — or as instructed by your dental professional.

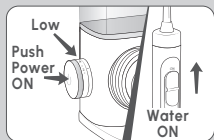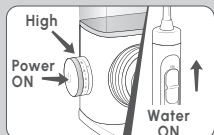Using the Water Flosser

Always place the tip in your mouth before turning the water ON. With water switch in OFF position, turn power on. **Lean over sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink.** Turn water switch ON. **When finished, turn water switch and power OFF.** During use you can stop the flow of water by sliding the water switch on handle to the OFF position.

NOTE: Motor will continue to run.

Recommended Technique

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

When Finished

Turn the unit OFF. Empty any liquid left in the reservoir. Unplug the unit.

Note: If you used mouthwash or other solution besides water, see "Using Mouthwash and Other Solutions" section.


NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE


Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

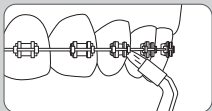Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

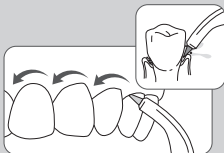Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.


To use the Pik Pocket™ tip, set the unit to the lowest pressure setting. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

CLEANING AND TROUBLESHOOTING YOUR WATERPIK™ WATER FLOSSER

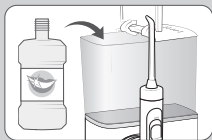Using Mouthwash and Other Solutions

Your Waterpik™ Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty. **Never immerse in water or any other liquid.**

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior. It is recommended that the reservoir is removed and cleaned weekly. The water reservoir is removable for easy cleaning and is top rack dishwasher safe heat cycle off. **Using toothpastes containing peroxide may damage the plastic housing of this product.** Using abrasive or chemical cleaners may damage the plastic housing of this product.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Run unit until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance. If using mouthwash or other solutions, more frequent cleaning is recommended.

Leakage

If leaking due to a damaged hose, replacement hoses and handles can be purchased at www.waterpik.com (US) or www.waterpik.ca (Canada) or contact Water Pik Inc. customer service for a self service repair kit.

Service Maintenance

Waterpik™ Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to www.waterpik.com (US) or www.waterpik.ca (Canada). You may also call our toll-free customer service number at 1-800-525-2774 (US) or 1-888-226-3042 (Canada). Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Still have questions? Go to www.waterpik.com (US) or www.waterpik.ca (Canada).

BATTERY CARE AND BATTERY DISPOSAL/RECYCLING

Battery Care

If the Toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage.

Do not expose the battery or appliance to excessive temperatures.

Battery Disposal/Recycling





This product contains a **rechargeable lithium ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

Please contact your local waste administrative department for any questions.

Note: The rechargeable battery is not replaceable. Remove only at end of product life. This process is not reversible and voids the product warranty.

BATTERY REMOVAL

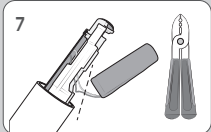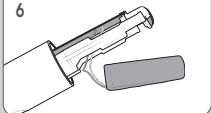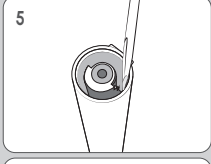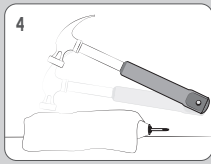**Battery Removal**

To remove the rechargeable battery for recycling, please discharge the battery fully before removing. Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Run down the battery by operating the Toothbrush Power Button until the motor does not run.
2. Obtain a towel or cloth, hammer, flat-head screwdriver, and wire cutters.
3. Remove and discard the brush head. Cover the handle entirely with the towel or cloth.
4. Clamp the top of the handle against your work surface with one hand and strike the handle housing around the bottom perimeter on all sides to release the end cap latch. Continue striking and rotating the handle until the end cap releases.
5. Remove the internal components by turning the handle upside down, then un-latch the plastic tabs and slide the chassis from the outer housing.
6. The internal components should slide out when the handle is positioned upside-down. If this does not occur, repeat steps 5 and 6 until the components release.
7. Remove the battery from the chassis. Cut the battery wires **one at a time** with the wire cutters.
8. Recycle and dispose of the battery properly. The remaining product pieces should be discarded appropriately.

Recycle/dispose of the battery at your nearest recycling center or take to a hazardous waste collection facility to comply with local and state regulations.

WARNING:

DO NOT cut both battery wires at the same time.

DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this product that it is free from defects in materials and workmanship for 2 years from the date of purchase. **Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim.** Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips and brush heads. **Use with improper voltage will damage product and VOID warranty.**

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

Arbitration

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik products or the general public. Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

FCC Statement (USA)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by Water Pik, Inc. could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by

turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

IC Statement (Canada)

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

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Church & Dwight Co., Inc.